


Hettich
BENELUX
LABORATORY EQUIPMENT



HYPOTHERMIA

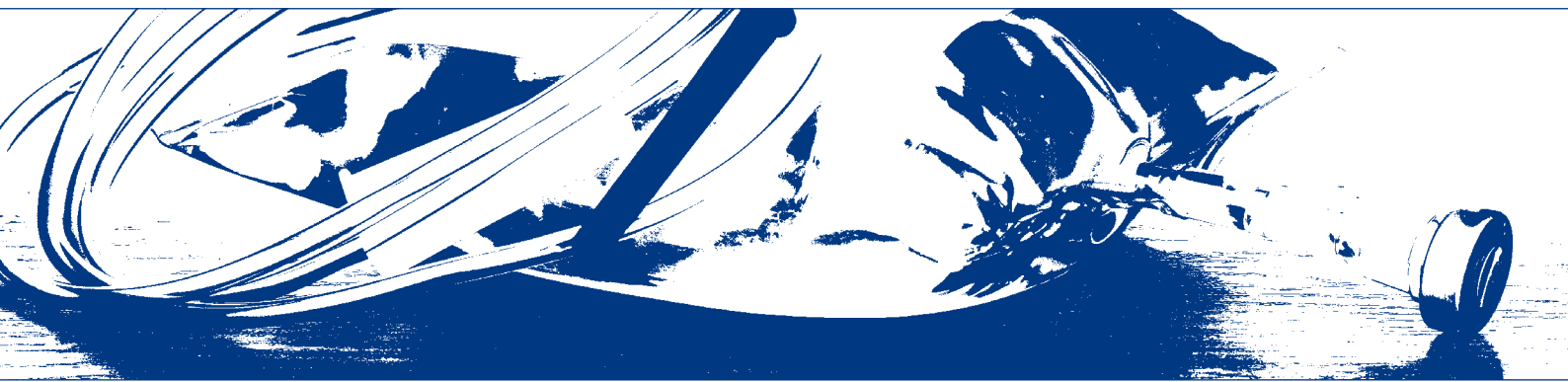

ENGINEERED
FOR YOU

Your leading supplier of climate control equipment



Hettich Benelux is your leading supplier of climate control equipment, which is tailored to your specific applications and requirements. We realise

that standard solutions are not enough in many cases, which is why Hettich Benelux continues to develop individual solutions for specific situations. Engineered for you.



Quality, innovation and service

Expect quality, innovation and service from Hettich Benelux, in everything we do. You will see it in our products, A-class parts, certifications, eco-awareness, service and communication. Not only is innovation important, the maintenance of your equipment is also in good hands.

Hettich Benelux is responsible for the design-to-production process; full, extensive tests are conducted on the end product. This applies to all incubators, ovens, refrigerators and freezers. Our products do not only meet the current requirements, but also anticipate future demands.



Our Hypothermia chambers are used worldwide, for example in hospitals, operating theatres, and Accident and Emergency departments (A&E).



Engineered for you

The hypothermia/heating chamber is used for storing intravenous fluids, dialysis fluids, irrigation fluids, filled syringes, (synthetic) blankets, hand towels and medication. We can also modify the equipment to

suit your specific requirements. As we develop and produce the hypothermia/heating chamber ourselves, we are able to implement the required modifications in a simple and practical way.



Flexible arrangement of chamber

Chambers can be customised to be built into the wall or free-standing chambers can be made with single or double doors, made dustproof from floor to ceiling. Combination chambers are also available,

with compartments that work independently of each other. Any combination of heated and chilled compartments can be manufactured. The hypothermia/heating chambers are designed so the user can choose shelves made of wire

mesh, sealed stainless steel shelves, and drawers of different depths, with or without compartments. Any hypothermia/heating chamber can be adapted to meet your requirements and area of application.



The benefits of Hettich Benelux

- Can be supplied as free-standing or built-in customised chamber with drawers or doors.
- Low energy consumption.
- Quick recovery times.
- The chambers are easy to open thanks to their ergonomic design.
- Visual alarm, easily visible due to red and green LED lighting.
- Easy to clean due to fully stainless steel interior.
- Adjustable interior arrangement.
- Self-closing door and drawers.
- Dust filter (optional) for use with blankets/flannels.



Hettich Benelux

Stable and precise temperature regulation | Clear visual alarm
environmentally friendly and energy efficient | Customised
Guaranteed quality

Energy-efficient with respect for the environment

The simple control system allows you to accurately set and read the temperature. The programmed 'neutral zone' ensures maximum energy-saving. PID sensors means a proportional cooling and heating until the set value is reached followed by an energy-neutral interval during

which time no cooling or heating takes place. This technique does not only guarantee a constant temperature, but it is also many times more eco-friendly than conventional methods.

Guaranteed quality

The Hettich Benelux quality management system meets the ISO 9001, ISO 14001


and ISO 13485 standards. The quality of the cabinets is guaranteed by measuring the temperature variation and fluctuation. A calibration certificate will be delivered with the equipment.

This product is registered with the Dutch Ministry of Health, Welfare and Sport (VWS) as a Class I medical appliance under name and number: Hettich ESP Hetherm (NL-CA002-2012-25728)



Engineered for you

Your Hettich Benelux partner

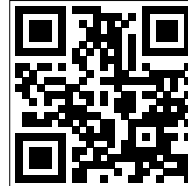

Member of SEDEX

Hettich Benelux B.V. | De Aaldor 9 | 4191 PC Geldermalsen | The Netherlands
www.hettichbenelux.com | +31 (0)88 221 99 00 | info@hettichbenelux.com

Vraag naar de specificatiefolders.  @HettichBenelux  Hettich Benelux